F.H.P.

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elementary English version
important content translated

Standing frames and rehabilitation seats
all presented products are intended for use with kids

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
 „ELEPHANT” STANDING FRAME is designed for an active rehabilitation of disabled children. The standing frame allows for independent (but under supervision) and safe standing in the appliance, fun and education. Thanks to its simplicity, the „ELEPHANT” STANDING FRAME is an extremely precious device at the rehabilitation of disabled children at home, in hospitals and rehabilitation centers. The correct use of this appliance can help to improve the quality of children's life.

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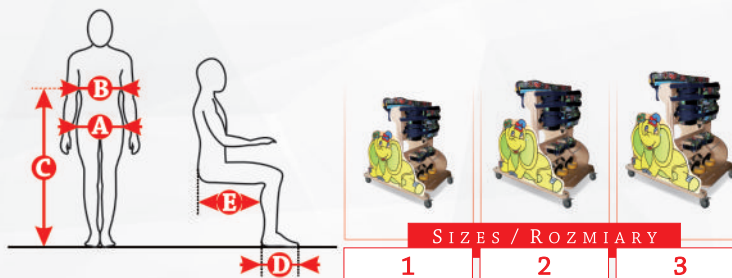ADJUSTMENT OF HIP CUSHION

In order to adjust the width of cushions supporting the chest, you need to loosen the knob, move the cushion to the desired position and then tighten the knob.

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MOVING

The device is equipped with 4 pieces of wheelset. Wheelsets have a ratchet brake.

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CHEST BELT AND HIP BELT

Thanks to the use of strips, belts are easy to fasten.

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CUSHIONS SUPPORTING CHEST AND HIPS

In order to regulate the width of cushions supporting the chest, you need to loosen knobs (located behind the standing frame), adjust cushions to the child's chest and tighten knobs. The position of hip cushions is similarly adjusted.

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SAFETY BELTS

Adjustment of the safety belt's length is realized by moving the tape in the buckle.

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ADJUSTMENT OF TABLE'S INCLINATION ANGLE

By loosening the nuts you can change the table's inclination angle. After setting a desired position, you need to tighten the nuts.

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VARIATION OF THE STANDING FRAME

In the electric version, the standing frame is adjusted by pressing an appropriate button on the remote control, and in the manual version by turning a proper knob, setting the standing frame and tightening screws.

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LOWER SUPPORT

This support is used to give the child a full support.

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
FOOTERS

The position of footers on the footrest plate is adjusted with the use of knobs.

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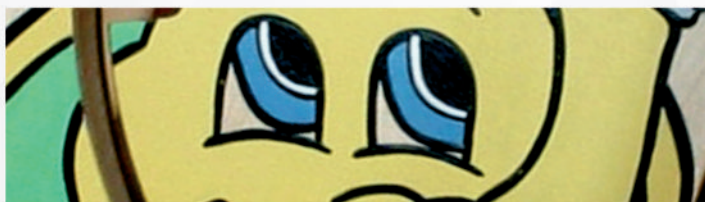
„ELEPHANT” STANDING FRAME is designed for an active rehabilitation of disabled children. The standing frame allows for independent (but under supervision) and safe standing in the appliance, fun and education. Thanks to its simplicity, the „ELEPHANT” STANDING FRAME is an extremely precious device at the rehabilitation of disabled children at home, in hospitals and rehabilitation centers. The correct use of this appliance can help to improve the quality of children’ life.

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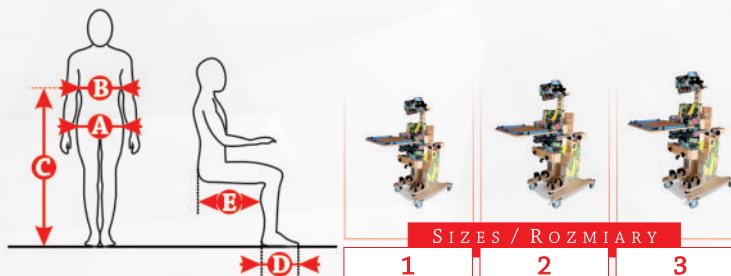SIZES / ROZMIARY			
	1	2	3
A. Width of hips [cm]	15 - 22	17 - 28	27 - 38
B. Width of the chest [cm]	15 - 22	17 - 28	27 - 38
C. Height from feet to armpit [cm]	48 - 72	58 - 88	68 - 98
D. Foot length [cm][max]	17	20	24
E. Distance from the back to the calf [cm]	11 - 23	19 - 35	24 - 43
1. User weight [kg][max]	25	45	60
2. Dimension of the device [cm][LxWxH]	80x60x65	86x65x75	92x66x80
3. Machine weight [kg]	25	30	35

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„DOG” rehabilitation standing frame is produced for children with disorders preventing the independent walking or sitting. It is used to the initial verticalization of children in the sitting, standing or lying position. Thanks to its multi-functionality, rich regulation and accessories, this device complexly completes the young patient's rehabilitation process. The standing frame is made entirely from eco-friendly and environmentally friendly materials. The use of wood, colorful upholstery and cheerful figure of a little dog in this frame causes that children more willingly use this appliance. The standing frame is built from a massive construction, and simultaneously it is very easy to use. The construction protects the child very well in the sitting, standing and lying position. This device can be used for fun and is a great addition in the rehabilitation process. Thanks to the use of chassis, the seat is mobile what allows the movement of the child, and in the case of tiredness – the child can rest after the rehabilitation in the lying position.

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SEAT INCLINATION ANGLE

It is possible to change the seat inclination angle, and with the complete disposition of the backrest and adjustment of the seat's angle, you can comfortably put the child in a lying and standing position.

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DEPTH OF THE SEAT AND BACKREST ANGLE

Children can sit in the seat with straight legs. They can also rest on the lying position after tiredness caused by the rehabilitation and fun.

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ABDUCTION BLOCK

If a rehabilitator decides that it is necessary to dismantle the abduction block, we can easily remove it.

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MOVING

The seat has a wheelset with a ratchet brake and steering wheels, so you can transport the child to another room.

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ADJUSTMENT OF FOOTRESTS

Each footrest can be adjusted separately depending on the child's needs and it is possible to change the inclination angle and height of footrests.

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SIDE PELOTES SUPPORTING CHEST AND HIPS

Adjustable side supports support the chest and hips, and simultaneously help to correct postural problems.

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TABLE

The table is included in the standard equipment of the „DOG” seat. It is possible to adjust the height and inclination angle of his device.

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ADJUSTMENT OF THE HEADREST

Adjustable headrest gives the possibility to set the height and width. There is also head suspender belt, which increases the safety and comfort of the child.

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KNEE SEPARATOR

Knee separator helps in the support of child's knees and their stabilization.

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ADJUSTMENT OF ARMRESTS

Inclination of the armrest can be adjusted in two ways:

- 1) Adjustment of the armrest, Adjustment of armrest's guide
- 2) Thanks to these adjustments, we can raise the young patient's comfort.

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
„SUN” STANDING FRAME is designed for an active rehabilitation of disabled children. The standing frame allows for independent (but under supervision) and safe standing in the appliance, fun and education.

Thanks to its simplicity, the „SUN” STANDING FRAME is an extremely precious device at the rehabilitation of disabled children at home, in hospitals and rehabilitation centers. The correct use of this appliance can help to improve the quality of children’ life.

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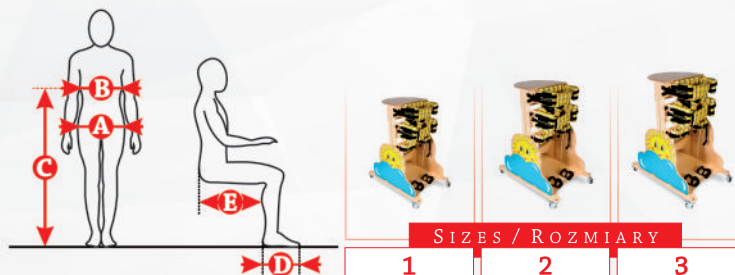	1	2	3
A. Width of hips [cm]	15 - 22	17 - 28	27 - 38
B. Width of the chest [cm]	15 - 22	17 - 28	27 - 38
C. Height from feet to armpit [cm]	48 - 72	58 - 88	68 - 98
D. Foot length [cm][max]	17	20	24
E. Distance from the back to the calf [cm]	11 - 23	19 - 35	24 - 43
1. User weight [kg][max]	25	45	60
2. Dimension of the device [cm][LxWxH]	80x60x65	86x65x75	92x66x80
3. Machine weight [kg]	25	30	35

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

„TURTLE and ELECTRIC TURTLE” rehabilitation standing frame is produced for children with disorders concerning the impossibility of independent sitting or walking. It is used to the initial verticalization of children in sitting, standing and lying position. It completes the rehabilitation process of the young patient via its multifunctionality, rich regulation and accessories. The standing frame is entirely made from eco-friendly and environmentally friendly materials. The use of wood, colorful upholstery and a cheerful figure of TURTLE causes that children more willingly use this appliance. The standing frame is built from a massive construction, and simultaneously it is very easy to use. The construction protects the child very well in the sitting, standing and lying position. This device can be used for fun and is a great addition in the rehabilitation process. Thanks to the use of chassis, the seat is mobile what allows the movement of the child, and in the case of tiredness – the child can rest after the rehabilitation in the lying position.

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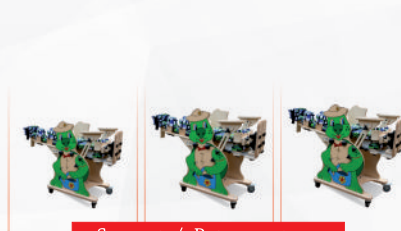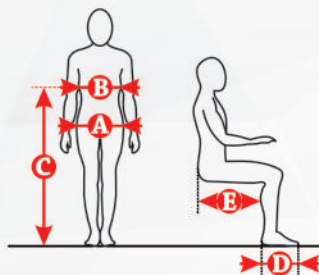	1	2	3
A. Width of hips [cm]	15 - 22	17 - 28	27 - 38
B. Width of the chest [cm]	15 - 22	17 - 28	27 - 38
C. Height from feet to armpit [cm]	48 - 72	58 - 88	68 - 98
D. Foot length [cm][max]	17	20	24
E. Distance from the back to the calf [cm]	11 - 23	19 - 35	24 - 43
1. User weight [kg][max]	25	45	60
2. Dimension of the device [cm][LxWxH]	80x60x65	86x65x75	92x66x80
3. Machine weight [kg]	25	30	35

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
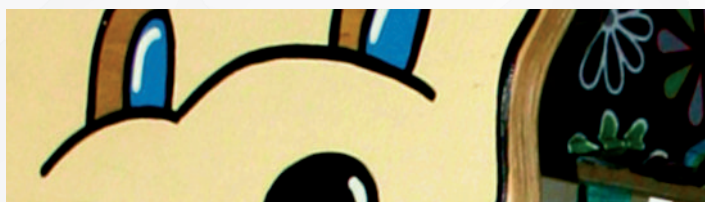
„TEDDY BEAR” seat is designed for an active rehabilitation of disabled children. This seat causes that children more quickly explore the world around them. „TEDDY BEAR” seat is used for the initial verticalization in sitting and lying position. Moreover, it ensures the full stabilization in both positions. It is perfect for both spastic children and children with reduced muscle tension. „TEDDY BEAR” seat allows for playing and education. It is a device that can be used both in hospital rehabilitation centers and at home. „TEDDY BEAR” seat is very useful in orthopedic disorders. The seat has a wheel steers allowing the movement of the child.

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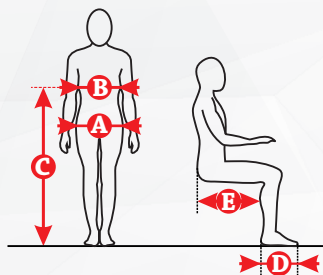	1	2	3
A. Width of hips [cm]	15 - 22	17 - 28	27 - 38
B. Width of the chest [cm]	15 - 22	17 - 28	27 - 38
C. Height from feet to armpit [cm]	48 - 72	58 - 88	68 - 98
D. Foot length [cm][max]	17	20	24
E. Distance from the back to the calf [cm]	11 - 23	19 - 35	24 - 43
1. User weight [kg][max]	25	45	60
2. Dimension of the device [cm][LxWxH]	80x60x65	86x65x75	92x66x80
3. Machine weight [kg]	25	30	35

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
Rehabilitation seat "SEAL" is irreplaceable in the initial verticalization. Thanks to the simple construction, it is easy to use, and simultaneously it perfectly fulfills its function. It can be safely used by very flaccid children with a reduced muscle tension and spastic children. Moreover, the seat can be very helpful in rehabilitation of orthopedic disorders. Thanks to the regulation of side pelottes, it helps to correct the postural problems. We give you the rehabilitation seat "SEAL", which is an extremely useful device in the rehabilitation of disabled children in home, hospital and rehabilitation center conditions. The proper use of this appliance can contribute to a significant improvement of children's life and facilitate their fun and education.

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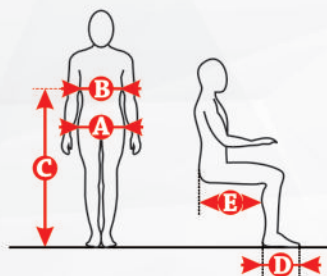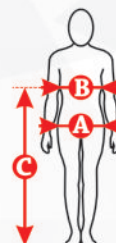
	1	2
A. Width of hips [cm]	15	17
B. Width of the chest [cm]	15	17
C. Height from feet to armpit [cm]	8	10
1. User weight [kg][max]	25	45
2. Dimension of the device [cm][HxLxW]	50x23x52	55x27x57
3. Machine weight [kg]	5	7

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SIDE PELOTTES SUPPORTING CHEST AND HIPS

Adjustable side pelottes support the chest and hips, and simultaneously help to correct postural disorders.

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REGULATION OF HEADREST

Adjustable headrest gives the ability to set the height and width. Moreover, it is possible to add the belt supporting the head, which increases the child's comfort.

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TABLE

It is possible to adjust the height and the inclination angle of the table what increases the comfort of a patient, who uses the rehabilitation seat "SEAL".

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SAFETY VEST

Safety vest has stripes. The length of these stripes is adjustable.

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REGULATION OF FOOTRESTS

It is possible to adjust footrests depending on the child's needs. The seat also has straps, which stabilize feet.

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SEAT DEPTH

Regulation of the seat depth depending on the child's needs.

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ABDUCTION BLOCK

Abduction block is an element, which can be easily removed if it is necessary.

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