
**i** Opis ogólny / Basic description

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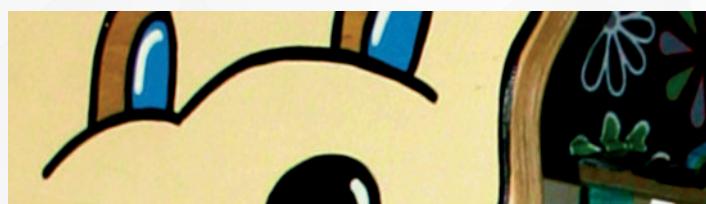
„TEDDY BEAR” seat is designed for an active rehabilitation of disabled children. This seat causes that children more quickly explore the world around them. „TEDDY BEAR” seat is used for the initial verticalization in sitting and lying position. Moreover, it ensures the full stabilization in both positions. It is perfect for both spastic children and children with reduced muscle tension. „TEDDY BEAR” seat allows for playing and education. It is a device that can be used both in hospital rehabilitation centers and at home. „TEDDY BEAR” seat is very useful in orthopedic disorders. The seat has a wheel steers allowing the movement of the child.

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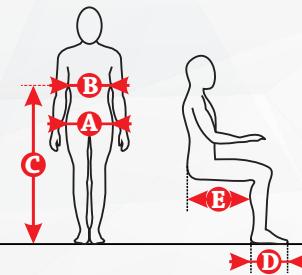## SEAT INCLINATION ANGLE

It is possible to change the seat inclination angle, and with the complete disposition of the backrest and adjustment of the seat's angle, you can comfortably put the child in a lying and standing position.

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## SIDE PELOTTES SUPPORTING CHEST AND HIPS


Ajustable side supports support the chest and hips, and simultaneously help to correct postural problems.

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## DEPTH OF THE SEAT AND BACKREST ANGLE

Children can sit in the seat with straight legs. They can also rest in the lying position after tiredness caused by the rehabilitation and fun.

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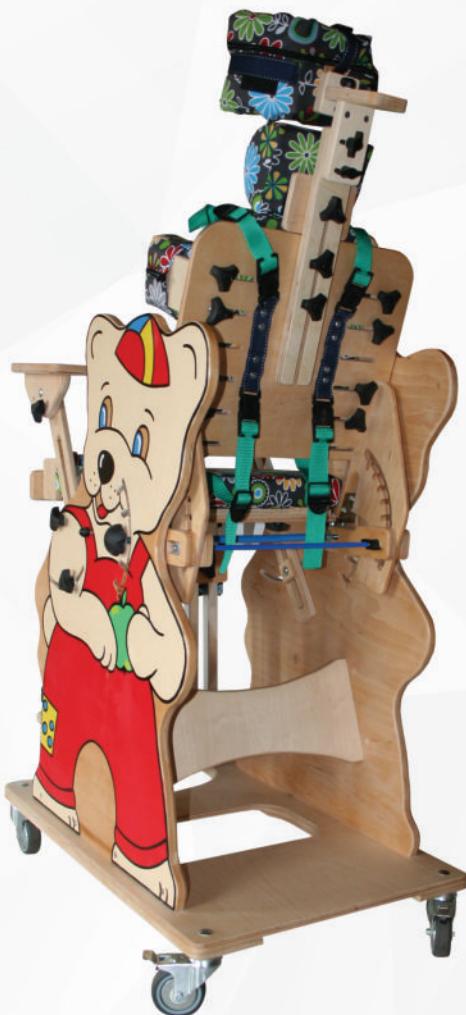## ADJUSTMENT OF ARMRESTS

The table is included in the standard equipment of the „DOG” seat. It is possible to adjust the height and inclination angle of his device.

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## MOVING

The seat has a wheelset with a ratchet brake and steering wheels, so you can transport the child to another room.

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## ADJUSTMENT OF THE HEADREST

Adjustable headrest gives the possibility to set the height and width. There is also head suspender belt, which increases the safety and comfort of the child.

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## ADJUSTMENT OF FOOTRESTS

Each footrest can be adjusted separately depending on the child's needs and it is possible to change the inclination angle and height of footrests.